

Dual Protection  
Safety I.V. Catheter

# Supercath™ 5

*A New Generation of  
Safety I.V. Catheter*



**MEDIKIT**


**Made in Japan**



# Needlestick Injuries ◀ **DUAL PROTECTION** ▶ Blood Exposures

## Safety Technology

### Needlestick Injuries & Blood Exposures

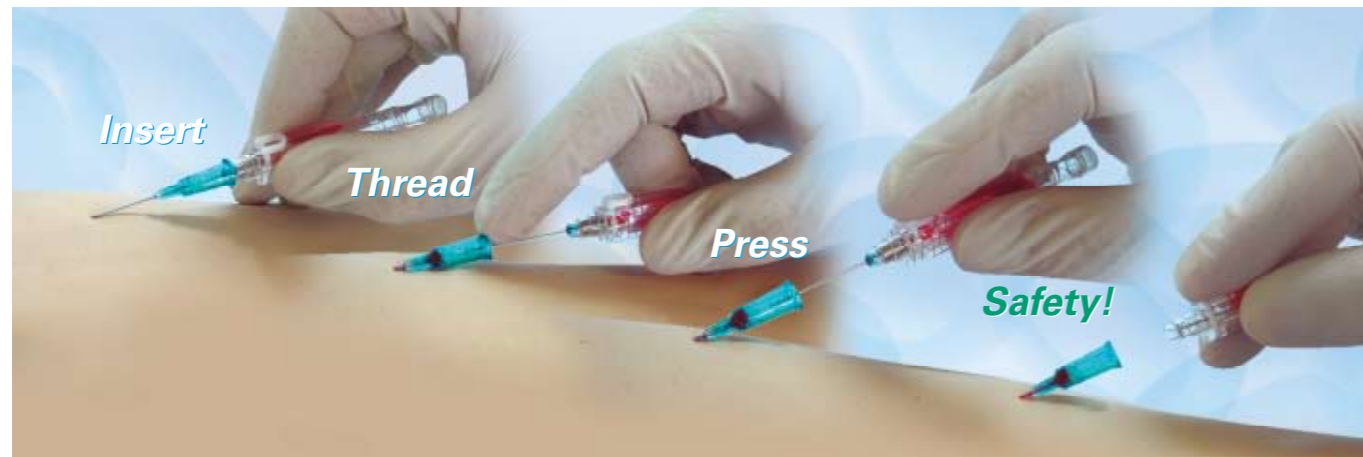



#### Fully Encased Needle Protection

- Reduces risk of accidental needlestick injuries to healthcare workers.
- Protects from blood exposure to contaminated needle.
- Offers a reassuring and confident feeling of safety.



#### Easy One-Handed Safety Activation

- When the safety mechanism is triggered, introducer needle is quickly retracted into needle shield.
- Minimizes in-service training.

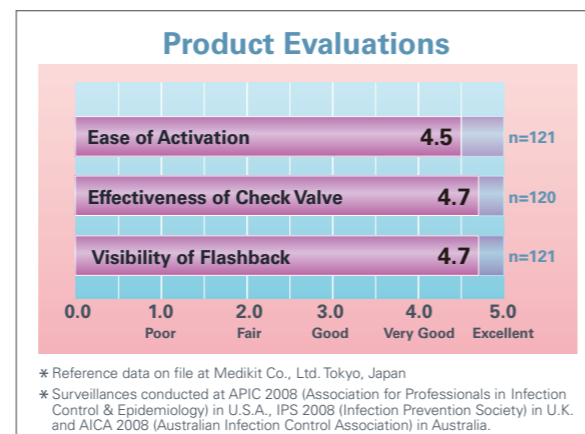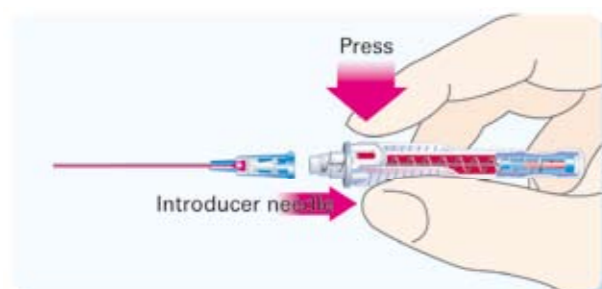

## Valve Technology

### Blood Exposures




#### Patented Integral Check Valve

- Designed to prevent blood exposure from catheter hub during venipuncture.
- Assists digital pressure to prevent blood spillage after withdrawal of the needle.

#### Non-Tortuous, Straight-Through Fluid Path

- Ensures a smooth fluid flow.

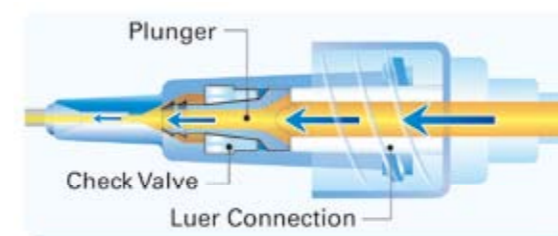

#### Catheter Hub Design Same as Conventional I.V. Catheters

- Simpler taping and patient comfort.


#### Latex Free & DEHP Free






- Eliminates dermal reaction to latex and eliminates a substance potentially harmful to patients.

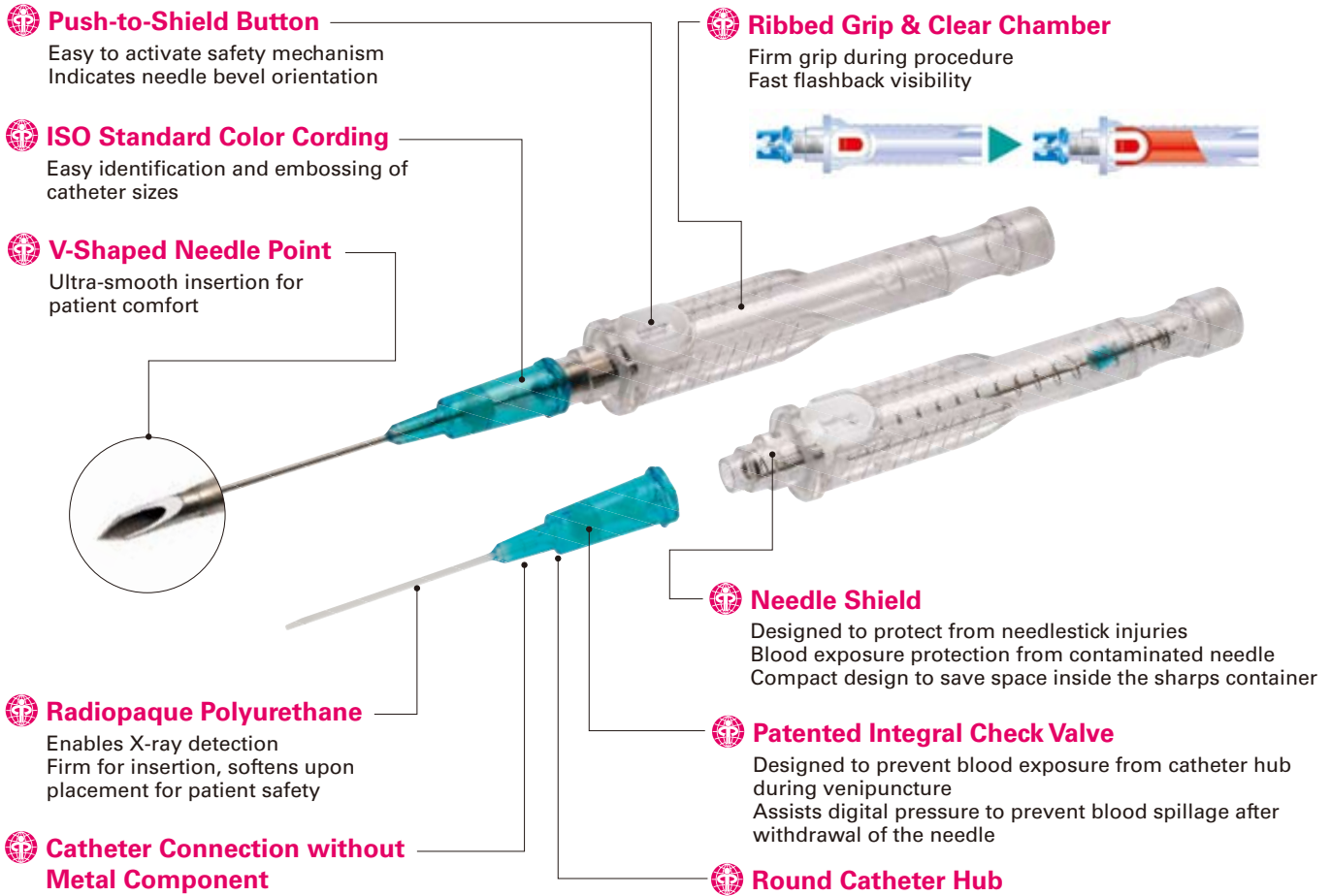

Dual Protection  
Safety I.V. Catheter

# Supercath™ 5



Color Code ISO10555-5:1996	Catalog Number	Gauge	Catheter			Flow Rate mL/min
			O.D.	I.D.	Length	
 Deep Green	SP120-18-31	18G	1.3 (mm)	0.9 (mm)	1 1/4" ( 31mm )	100
 Pink	SP120-20-31	20G	1.1 (mm)	0.7 (mm)	1 1/4" ( 31mm )	55
 Deep Blue	SP120-22-25	22G	0.9 (mm)	0.6 (mm)	1 " ( 25mm )	31
 Deep Blue	SP120-22-31	22G	0.9 (mm)	0.6 (mm)	1 1/4" ( 31mm )	28
 Yellow	SP120-24-19	24G	0.7 (mm)	0.4 (mm)	3/4" ( 19mm )	16

※ Please refer to the package insert for complete instructions for use.

50pcs/box, 20boxes/carton



**MEDIKIT Co.,Ltd.**

13-2, Yushima 1-chome, Bunkyo-ku, Tokyo, 113-0034 Japan  
Phone: 81-3-3839-8633 Fax: 81-3-3839-3978

[www.medikit.co.jp/english/](http://www.medikit.co.jp/english/)

Manufacturer

**TOGO MEDIKIT Co.,Ltd.**

CE 0197